



Support Surfaces

PM211AA

8" Comfortable Pressure Relief Mattress Replacement System

The Auto Adjust Alternating Pressure System is a mattress replacement designed with the very latest in microprocessor technology, allowing for simple operations while providing maximum patient care. This system automatically and accurately self-adjusts, providing optimal patient comfort based on the individual's weight. It relieves caregivers from manual configuration of the mattress.






Features:

- 3 modes of operation: Alternating Pressure, Static and Seat Inflation.
- Maximum Inflation Mode promotes a safe environment for patient transport and also facilitates medical intervention by medical staff.
- Control panel has a soft-touch keypad that is easy to wipe down and clean.
- Microprocessor control is ideal for automating procedures in order to simplify operations. It is also very durable and reliably consistent.



Mattress Cover has a full peripheral zipper and is made from an anti-bacterial, waterproof and vapor permeable material that reduces shear and friction.

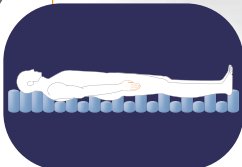

Alternating Pressure mode relieves pressure points beneath patient, preventing and managing pressure ulcers.



Convenient transport-mode cap maintains mattress air pressure during transport and prolonged power outages.



Visual/Audible alarms alert caregivers to low pressure or power failure conditions.



Lightweight and compact pump is quiet, energy efficient and ideal for transport and handling.

Specifications:

PUMP	Electrical	AC120V 60Hz 0.17A, AC230V 50/60Hz 0.07A	
	Air Output (lpm)	8 liter/min	
	Visual low pressure alarm	YES	
	Apex Quick Connector	YES	
	Static Function	YES	
	Cycle Time (min)	10 min	
	Dimension (L x W x H)	29.1 x 20 x 11.7 cm / 11.4" x 7.9" x 4.6"	
MATTRESS	Weight	2.2 Kg / 5 lb	
	Arrangement	20x8" Air Bladders	
	Material	Cover	Nylon/PU
		Cells	Nylon/PU
		Base	Nylon/PU
	Dimension (L x W x H)		200 x 90 x 20.3 cm / 78.7" x 35.4" x 8"
			200 x 85 x 20.3 cm / 78.7" x 33.5" x 8"
		200 x 80 x 20.3 cm / 78.7" x 31.5" x 8"	
Weight Capacity	180 kg or 400 lb		



Working Together With Care...

Distributed by:

IHS
SALES & SERVICE

Independent Healthcare Service Ltd.
Unit 3 Westpoint Business Park,
Mullisduart, Dublin 15.

Tel: 01 443 3390
Email: info@ihsltd.ie
Website: www.ihsltd.ie

Working Together With Care...