



RELAX - INNOVATIVE PERFECTION

The GK-RELAX side entry bath is designed to be used in many different bathing areas. An electrically powered wide opening door allows quick and easy side access.

The design of the bath also allows easy access for most patient hoists.



GK Bath Tubs Technology
RELAX.



	SALES & SERVICE
	Tel: 01 443 3390
	Email: info@ihsltd.ie
	Website: www.ihsltd.ie
<i>Working Together With Care...</i>	

Independent Healthcare Service Ltd.
Unit 3 Whospine Business Park,
Malahide, Dublin 15.

REINDEL
LANKER

Typing, spelling or transmission errors accepted for all contents. All product related information listed and binding and is subject to change for technical production or other reasons (date: 01/2020).



RECLINING-TUB RELAX






6 Filling of the bath, adjustment of disinfectant and hydro-massage are all controlled from the control panel. The water temperature is controlled by an approved thermostatic valve.

7 Fill in spout with integrable inlet of bath additive.

8 The hydro-massage system is a fully adjustable patented 2 jet system. It is self emptying and can be disinfected.

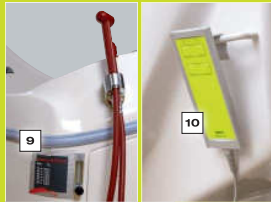

9 The GK-Relax is available with an electronically controlled disinfection system. This allows the amount of disinfectant to be adjusted to the required level.

10 Remote control



DETAILS.

- 1** Electric door locking
- 2** Grip bar
- 3** Drain assy
- 4** Swivel entrance door
- 5** Integrated head rest

The Relax Bath has an automatic fill, which allows a set volume of water at a controlled temperature to fill the bath. This saves time allowing more time for patient care.

TECHNICAL DATA:	
Supply:	<ul style="list-style-type: none"> • Hot and cold water with check valve 1/2" • Connector DN 50 • Power supply 230V/50Hz • equipotential bonding
Weight:	• ca. 150 kg
Tub filling:	• ca. 90 - 130 liters
Certificate:	• CE-Sign

