

CI/SIB	(74.28)	X	(U35)

Showering and Bathing



Shower Trolleys and Stretchers

The complete Showering & Bathing solution for the home, nursing home and hospital



Showering and Bathing

Shower Trolleys and Stretchers

The complete Showering and Bathing solution for the home, nursing home and hospital

The Health and Safety Executive has stated that handling heavy or awkward loads can result in musculoskeletal disorders, including back ache, pains, and strains.

The Shower Trolley range allows a patient to be transferred to the bathing area, showered, dried and clothed, minimising the need to lift manually and in line with current Moving and Handling Regulations.



T11 Shower/Changing Stretcher

The Shower/Changing Stretcher allows patients to be washed over a bath or in a shower area, in the prone or raised position. For ease of storage the stretcher folds flat against the wall when not in use.



T1 Shower Trolley


Heavier patients are easily accommodated, up to 190kgs (30 stones) and to aid transfer a simple foot pedal operation allows the height of the trolley to be adjusted to a convenient height.



T12 Stretcher

Shower Trolleys and Stretchers

+ T1 Shower Trolley shown

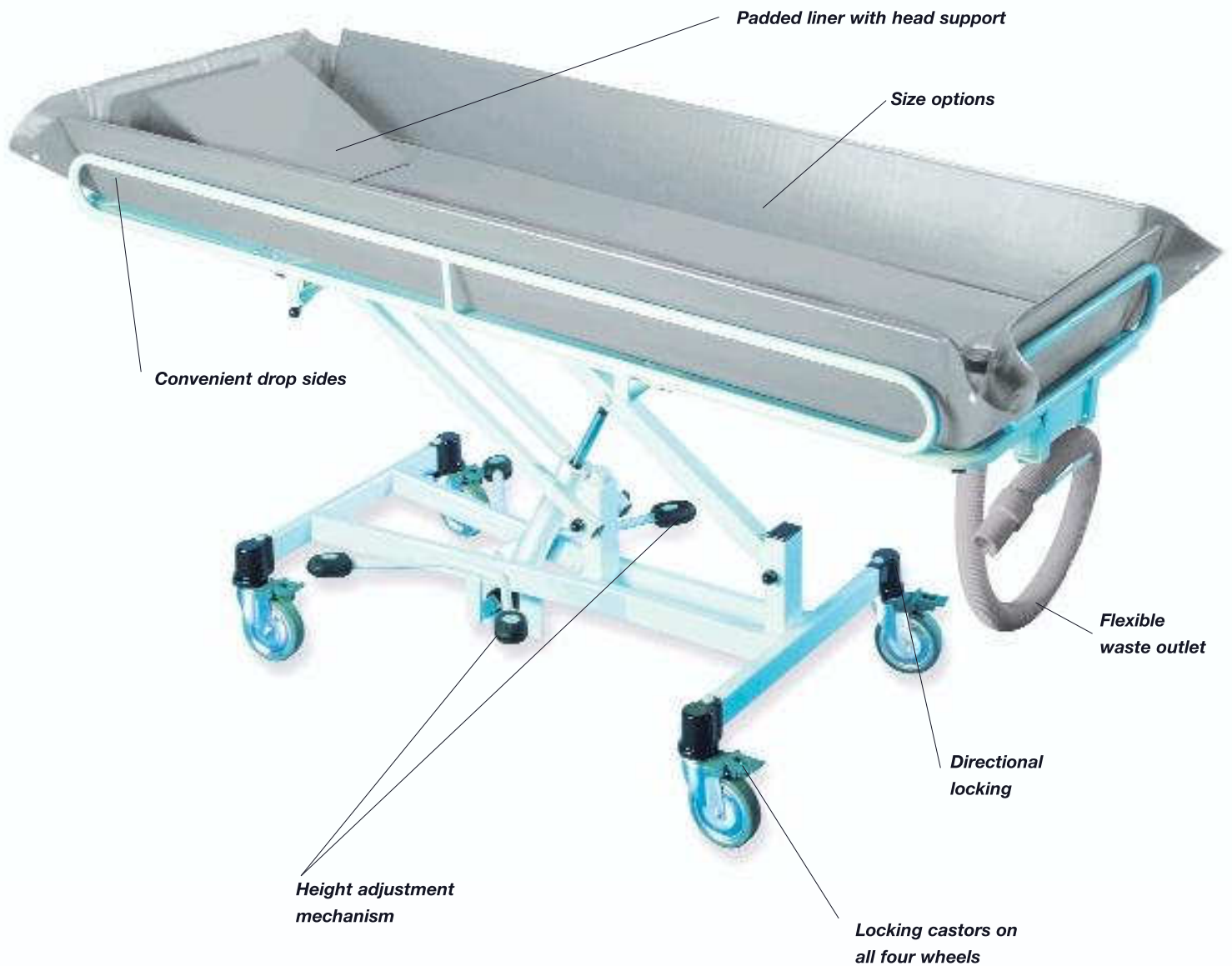

◆ T12 Shower/Stretcher shown



Shower Trolleys and Stretchers Features

Shower Trolleys

Easy steer mechanism

Makes moving between rooms or hospital wards easy.

Braked castors

Three locking castors prevent movement during patient transfer. The fourth castor has directional locking to assist manoeuvrability.

Drop sides

Convenient drop sides on both sides of the trolley to aid transfer.

Adjustable height

The height adjustment mechanism allows a carer to move the trolley up and down to a convenient height from any side, avoiding the need for stretching or leaning.

Comfort - Padded liner complete with head support for added comfort.

Reduced manual handling - Users can be showered and dried, minimising the need to lift manually.

Model options

+ **Full size model** - T1 Shower Trolley (1930mm x 750mm) will comfortably accommodate a large adult. The T2 is a fixed height trolley. The T4 is electrically operated by hand control

● **Smaller top model** - T3 Shower Trolley (1600mm x 700mm) with smaller top available where space is limited.

Shower/Changing Stretcher T11,12 & 13

The Chiltern Invadex Shower/Changing Stretcher range is safe and sturdy, providing reassurance for the bather. Choice of height helps to make showering easy for the carer.

Model options

- **T11** - Shower/ Changing stretcher.
- ◆ **T12** - with adjustable back can be locked in various positions to suit the user.
- **T13** - Shower/Changing unit features a padded liner with head support and drop sides.




T13 Shower/Changing unit
Fold away wall mounted model



- **Easy to store** - Folds neatly against the wall when not in use
- **Easy drain** - Drain holes allow water to run away quickly.





Shower Trolleys and Stretchers

Reducing the need for manual handling

The Chiltern Invadex range of Shower Trolleys can make a valuable contribution towards a safe moving and handling policy, helping to avoid back injury. The correct equipment, coupled with appropriate training can minimise the need for manual handling. Please ask for a copy of our comprehensive range of Moving and Handling training programmes.

Safety

Safe, secure transfer is ensured by a host of safety features. All models are manufactured to ISO 9001 quality standards.

Survey and Installation

A Chiltern Invadex Technical Sales Representative can provide a detailed survey. Please ask for details.



Chiltern Invadex Ltd
 Chiltern House
 6 Wedgwood Road
 Bicester, Oxfordshire
 OX26 4UL
 Tel: 01869 246470
 Fax: 01869 247214
 E-mail: sales@chilterninvadex.co.uk
 Web: www.chilterninvadex.co.uk

Overseas Enquiries
 Tel: +44 (0) 1869 365545
 Fax: +44 (0) 1869 365567
 E-mail: export@chilterninvadex.co.uk

SALES & SERVICE

IHS

Independent Healthcare Service Ltd.

Tel: 01 443 3390

Email: info@ihsltd.ie

Unit 3 Westpoint Business Park,
Mulhuddart, Dublin 15.

Website: www.ihsltd.ie

Working Together With Care...

Chiltern Invadex Ltd
 Unit 16 Evans Business Centre
 68-74 Queen Elizabeth Avenue
 Hillington Industrial Estate
 Glasgow
 G52 4NQ
 Tel: 0141 880 9986
 Fax: 0141 880 9990



Bicester is BS EN ISO 9001 approved (FM 40009).
 Bangor Northern Ireland is BS EN ISO 9002 approved (FM 21085).

Copyright © Chiltern Invadex Ltd
 The design of the product shown on this brochure is vested in Chiltern Invadex Ltd and should not be copied or reproduced in anyway whatsoever without the express permission in writing of Chiltern Invadex Ltd.
 As part of our commitment to new product development we reserve the right to alter specifications without prior notice.
 This is an original design of Chiltern Invadex Ltd in which Design Right subsists.
 Published by: Chiltern Invadex Ltd, 6 Wedgwood Road, Bicester, Oxon, OX26 4UL Registered in England No: 1182024 ISSUE 5 July 2003
 In Northern Ireland: Chiltern Invadex (NI) Ltd, 5 Balloo Way, Balloo Industrial Estate, Bangor, County Down, BT19 7QZ Registered No: NI 19958

Shower/Changing Stretcher
 Wall Mounted tip-up models

T12 Wall mounted tip-up with adjustable back

Plan view T12/T11

T11 Wall mounted tip-up

Shower Trolley
 Mobile and Wall mounted tip-up models

T1 Wheeled shower trolley

Plan view T1/T13

T13 Wall mounted tip-up shower bathing/changing unit

Technical Information

	A	B	C	D	E	F	G	H	J	K
T1 Trolley	-	700	-	-	1930	200	580	960	-	-
T2 Trolley	-	700	-	-	1930	200	Fixed	960	-	-
T3 Trolley	-	700	-	-	1600	200	580	960	-	-
T4 Trolley	-	750	-	-	1930	200	580	960	-	-
T11 Stretcher	132	700	-	900	1900	-	-	-	700	110
T12 Stretcher	132	700	ADJ	900	1900	-	-	-	700	110
T13 Stretcher	214	715	-	900	1910	200	-	-	750	185

Handing of T12 Shower Stretcher
 The adjustable backrest on the T12 is on the carers right hand side as standard.
 Opposite hand made to order.
 To arrange for Chiltern Invadex to fit your selected stretcher please contact your nearest Divisional Office.

Dimensions are in mm and for guidance only



Other products available from Chiltern Invadex include:

Showering and Bathing

- Shower Cubicles
- Level, Easy and Step-in Showers
- Shower Accessories
- Shower Seats
- Walk-in and Easy Transfer Baths
- Shower Chairs

Moving and Handling

- Ceiling Mounted Hoists
- Wheelchairs
- Slings
- Mobile Hoists
- Poolside Hoists

Brochure Hotline: 0800 917 0399