



Advance

Get closer



OXFORD

INNOVATE.
BUILD FOR THE FUTURE



Outstanding performance and lifting range

Advance

Oxford's design formula has taken into account the following key elements;

Ergonomic Considerations: making the device simple and safe to use.

Injury Prevention focus: ensuring that we add features that facilitate safe transfer practice

Human Factors: paying attention to the positioning of key interaction points such as handles and battery pack



Award winning Product

All of the above has delivered what we feel is truly and outstanding product concept. One reason why the Oxford Advance is the proud **winner** of a '**commendable design concept award**' from the Independent Living Foundation.



Advance

A compact product that gets you closer

The Advance's unique 'swan neck' leg design allows the hoist to get even closer to the widest of obstacles. This is especially useful when encountering large chairs, bulky commodes and wheel chairs.

The 'swan neck' design also ensures that the client feels safe and secure during the transfer.




Small yet incredibly capable

The Advance is compact yet able to perform an outstanding range of transfers. Whether from the floor, on to a nursing or domestic bed, the Advance copes with ease.

Key dimensions:

Min. lift range:	390 mm
Max. lift range:	1,690 mm
Max. leg clearance:	115 mm
Safe working load:	155 kgs
Folded compact dimensions: Height X 430 mm, Length X 1,180 mm, Width X 550 mm	

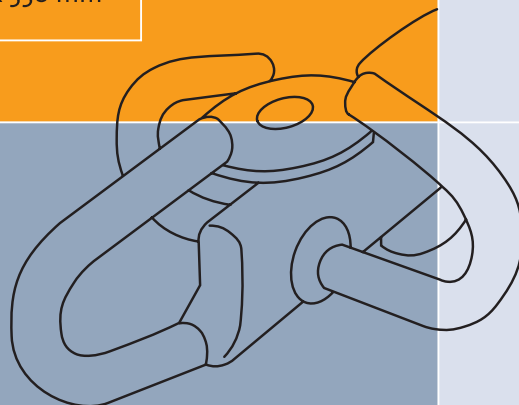

Exceptionally versatile

The Advance patient hoist from Oxford signifies a real stride forward in patient lift design and performance. It uses the very latest design and manufacturing technologies to ensure the Advance is the most functional and stylish hoist in its class today.

The Advance also has truly outstanding performance. The moment you move or operate the Advance you are assured that every little detail has been considered during its meticulous design.



With this level of performance you would expect the Advance to be large and cumbersome. However, this could not be further from the truth, the Advance is both small and nimble.

Furthermore, the Advance lends itself to those environments where space and storage are at a premium. It is especially useful in home settings and small nursing institutions, where space sometimes is simply not available.

Additional useful features



The oversized handle provides a large surface area from which the carer can position the lift before use.



Triangular folded design: This useful feature ensures two very important functions. (1) The hoist stands unaided when being stored or transported. (2) The hoist can be safely left in the knowledge that it will not become unstable whilst not in use.



The push footpad has been designed to reduce the force needed to initiate movement when manoeuvring the stationary lift in a forward direction.

**INNOVATE.
BUILD FOR THE FUTURE**



Easy to transport and store



Unlike other products on the market the Advance needs NO tools or special attachments to achieve its compact folded position. This functionality gives active users the option to travel outside the confines of their own home. Truly bringing a new dimension to mobile patient handling equipment.

14.9 kg



Base

14 kg



Mast and boom

In addition to outstanding performance, the Advance can also be folded, or separated for onward transportation.



Advance

Oxford Advance Technical Specifications

Electric Model

Specification	Imperial	Metric
Safe Working Load	341 lbs	155 kgs
Maximum Overall Length	51"	1,300 mm
Minimum Overall Length	49"	1,250 mm
Maximum Overall Height	73"	1,860 mm
Minimum Overall Height	53.5"	1,360 mm

Folded Dimensions		
Height	17.1"	450 mm
Depth	46.4"	1,180 mm
Width	21.6"	550 mm
Spreader Bar Max. Height	66.5"	1,690 mm
Spreader Bar Min. Height (usable)	15.3"	390 mm
Height at Maximum Reach	46"	1,170 mm
Reach at Maximum Height	25.6"	650 mm
Reach at Minimum Height	13.7"	350 mm
Maximum Reach*	32.8"	835 mm
Turning Radius	55.9"	1,420 mm
Legs Open - External Width	42.5"	1,080 mm
Legs Open - Internal Width	39.3"	1,000 mm
Legs Closed - External Width	26.3"	670 mm
Legs Closed - Internal Width	22.4"	570 mm
Overall Height of Legs	4.5"	115 mm
Ground Clearance	1.1"	30 mm
Front Twin Castors	3"	75 mm
Rear Braked Castors	3.9"	100 mm

*Reach = centre of spreader bar to the front of the mast

Weights		
Mast, Base & Boom Assembly	63.1 lbs	28.9 kgs
Power Pack	2.8 kg	2.8 kgs
N/A Total	69.9 lbs	31.7 kgs
Base Assembly (not incl. battery)	32.8 lbs	14.9 kgs
Mast & Boom (not incl. battery)	30.9 lbs	14.0 kgs

Specifications may vary due to manufacturing tolerances

Warranty

The Advance comes with a 2 year warranty as standard. Please refer to your user manual for further information.

With thanks to the Kidderminster Treatment Centre for the use of their facility during the photo shoot.



SALES & SERVICE

Tel: 01 443 3390

Independent Healthcare Service Ltd.
Unit 3 Westpoint Business Park,
Malahide, Dublin 15.

Email: info@ihsltd.ie

Website: www.ihsltd.ie

Working Together With Care...

